
**Dentistry — Brackets and tubes for use in
orthodontics**

Médecine bucco-dentaire — Consolés et tubes utilisés en orthodontie

STANDARDSISO.COM : Click to view the full PDF of ISO 27020:2010



PDF disclaimer

This PDF file may contain embedded typefaces. In accordance with Adobe's licensing policy, this file may be printed or viewed but shall not be edited unless the typefaces which are embedded are licensed to and installed on the computer performing the editing. In downloading this file, parties accept therein the responsibility of not infringing Adobe's licensing policy. The ISO Central Secretariat accepts no liability in this area.

Adobe is a trademark of Adobe Systems Incorporated.

Details of the software products used to create this PDF file can be found in the General Info relative to the file; the PDF-creation parameters were optimized for printing. Every care has been taken to ensure that the file is suitable for use by ISO member bodies. In the unlikely event that a problem relating to it is found, please inform the Central Secretariat at the address given below.

STANDARDSISO.COM : Click to view the full PDF of ISO 27020:2010



COPYRIGHT PROTECTED DOCUMENT

© ISO 2010

All rights reserved. Unless otherwise specified, no part of this publication may be reproduced or utilized in any form or by any means, electronic or mechanical, including photocopying and microfilm, without permission in writing from either ISO at the address below or ISO's member body in the country of the requester.

ISO copyright office
Case postale 56 • CH-1211 Geneva 20
Tel. + 41 22 749 01 11
Fax + 41 22 749 09 47
E-mail copyright@iso.org
Web www.iso.org

Published in Switzerland

Foreword

ISO (the International Organization for Standardization) is a worldwide federation of national standards bodies (ISO member bodies). The work of preparing International Standards is normally carried out through ISO technical committees. Each member body interested in a subject for which a technical committee has been established has the right to be represented on that committee. International organizations, governmental and non-governmental, in liaison with ISO, also take part in the work. ISO collaborates closely with the International Electrotechnical Commission (IEC) on all matters of electrotechnical standardization.

International Standards are drafted in accordance with the rules given in the ISO/IEC Directives, Part 2.

The main task of technical committees is to prepare International Standards. Draft International Standards adopted by the technical committees are circulated to the member bodies for voting. Publication as an International Standard requires approval by at least 75 % of the member bodies casting a vote.

Attention is drawn to the possibility that some of the elements of this document may be the subject of patent rights. ISO shall not be held responsible for identifying any or all such patent rights.

ISO 27020 was prepared by Technical Committee ISO/TC 106, *Dentistry*, Subcommittee SC 1, *Filling and restorative materials*.

STANDARDSISO.COM : Click to view the full PDF of ISO 27020:2010

Introduction

This International Standard has been developed as a result of the difficulty often encountered by clinicians to make meaningful comparisons between brackets and tubes using the information currently available from manufacturers and suppliers.

STANDARDSISO.COM : Click to view the full PDF of ISO 27020:2010

Dentistry — Brackets and tubes for use in orthodontics

1 Scope

This International Standard is applicable to brackets and tubes for use in fixed orthodontic appliances.

This International Standard gives details of methods to compare the functional dimensions of orthodontic brackets and tubes, the test methods by which they can be determined, as well as packaging and labelling information.

This International Standard does not specify specific qualitative and quantitative requirements for freedom from biological hazards; which are covered in ISO 10993-1 and ISO 7405.

2 Normative references

The following referenced documents are indispensable for the application of this document. For dated references, only the edition cited applies. For undated references, the latest edition of the referenced document (including any amendments) applies.

ISO 1942, *Dentistry — Vocabulary*

3 Terms and definitions

For the purposes of this document, the terms and definitions given in ISO 1942 and the following apply.

3.1

angle of torque

θ

occlusal-gingival angle formed between the intersection of the line perpendicular to the tangent to the tooth side surface of the base and the line bisecting the slot in the occlusal-gingival direction, when viewed along the mesial-distal long-axis of the slot

NOTE If the angle is oriented toward the occlusal (gingival), it is designated as positive (negative).

See Figure 1.

3.2

angulation

α

angle between the line perpendicular to the mesial-distal long-axis of the slot and the central occlusal-gingival axis of the bracket/tube, when viewed from the buccal/labial

NOTE The angulation is positive (negative) when the gingival part of the central occlusal-gingival axis is inclined toward the distal (mesial) relative to the line perpendicular to the mesial-distal long-axis of the slot.

See Figure 2.

3.3

archwire

wire that is formed approximately to the shape of the dental arch

3.4

auxiliary slot

additional opening in the bracket or tube for insertion of auxiliary elements

3.5

auxiliary slot dimensions

diameter/cross-sectional dimensions of the largest wire that passes through an auxiliary slot

3.6

band

structural annular component affixed to the outer circumference of the crown of a tooth and to which a bracket or tube can be attached

3.7

base

part of the bracket or the tube that is attached to the enamel or band

3.8

descriptor

code to identify the nominal slot height in thousandths of an inch without unit designation, in accordance with accepted orthodontic practice (e.g. 18 or 22)

3.9

bracket

tube

structural unit attached to a band or base that is capable of retaining an archwire

3.10 in-out

3.10.1

bracket in-out

l_i

distance between the floor of the slot and the tooth side surface of the base along the intersection of the plane perpendicular to mesial-distal long-axis of the slot in the centre of the bracket slot and the plane bisecting the slot in the occlusal-lingual direction, when viewed along the mesial-distal long-axis of the slot

See Figure 3 a).

3.10.2

tube in-out

l_i

distance between the floor of the slot and the tooth side surface of the base, at the mesial end of the tube and in the midplane of the slot in the occlusal-lingual direction; for tubes with a mesial chamfer, the in-out is the distance at the mesial end of the tube between that part of the floor of the slot in contact with the largest dimension of an archwire (round wire for a round tube and a rectangular wire for a rectangular tube) able to pass the entire tube, and the tooth side surface of the base in the midplane of the slot in the occlusal-lingual direction, when viewed from the lingual

See Figure 3 b).

3.11 rotational offset

δ

angle between a line parallel to the floor of the slot and the line connecting the points of intersection of the lines along the mesial and distal end-faces of the slot at the tooth side surface of the base, when viewed from the occlusal

NOTE When the distal (mesial) part of the floor of the slot is furthest from the tooth it is known as a distal (mesial) offset.

See Figure 4.

3.12 slot

mesial-distal oriented opening in a bracket or tube primarily to contain an archwire

3.13 slot depth

d

minimum buccal-lingual dimension of a rectangle where the projection of its buccal/labial side is tangent to the shorter side of the intended slot, when viewed along the mesial-distal axis of the slot

See Figure 5.

3.14 slot height

h

maximum occlusal-gingival dimension of a rectangle that is fully engaged in the slot, when viewed along the mesial-distal long-axis of the slot

See Figure 5.

3.15 slot length

l_s

smallest mesial-distal dimension of the slot between the mesial and distal ends of the slot

4 Requirements

4.1 Hazardous elements

For the purposes of this International Standard, cadmium, beryllium, lead and nickel are designated to be hazardous elements and the manufacturer shall state the concentrations as a mass fraction expressed as a percentage.

4.2 Measurement of dimensions

4.2.1 The following dimensions shall be recorded to the nearest 0,01 mm. When determined in accordance with Clause 5, the following dimensions of the product shall comply with the ranges stated by the manufacturer:

- a) in-out l_i
- b) slot depth d
- c) slot height h
- d) slot length l_s
- e) dimensions of each auxiliary slot.

4.2.2 The following angles shall be recorded to the nearest 1°. When determined according to Clause 5, the following angles of the product shall comply with the ranges stated by the manufacturer:

- a) angle of torque θ
- b) angulation α
- c) the rotational offset δ

5 Test methods

5.1 Sampling

Six specimens of a single product from one batch shall be procured for each test.

Measurements shall be made on each dimension listed in 4.2 of each sample.

5.2 Dimensions

5.2.1 Apparatus

Measurements shall be taken with callipers, micrometers, optical comparators, or other devices with an accuracy of 0,005 mm or 0,5°.

5.2.2 Measurement procedures

5.2.2.1 Angle of torque (see Figure 1)


5.2.2.1.1 When viewed from the labial/buccal, perform the following.

- a) Construct a line along the mesial-distal long-axis of the slot (X-axis).
- b) Construct plane 1 (see A–A in Figure 1) perpendicular to the X-axis in the centre of the slot (depicted by the symbol \perp A).


5.2.2.1.2 When viewed in plane 1, perform the following in plane 1.

- a) Construct line 2 bisecting the slot in the occlusal-gingival direction.
- b) At the point of intersection of the tooth side surface of the base and line 2 (point 3 in Figure 1), construct line 4 tangent to the tooth side surface of the base.
- c) Construct line 5 perpendicular to line 4 (depicted by the symbol \perp B) passing through point 3.
- d) Measure the angle of torque, θ , between line 5 and line 2.

A suggested method is to view along the mesial-distal long-axis of the slot (X-axis) of the bracket or tube and place a convex radius template in plane 1 of the bracket or tube, against the tooth side surface of the base, thus constructing the tooth side surface of the base. Then construct lines 2, 4 and 5 as defined above.

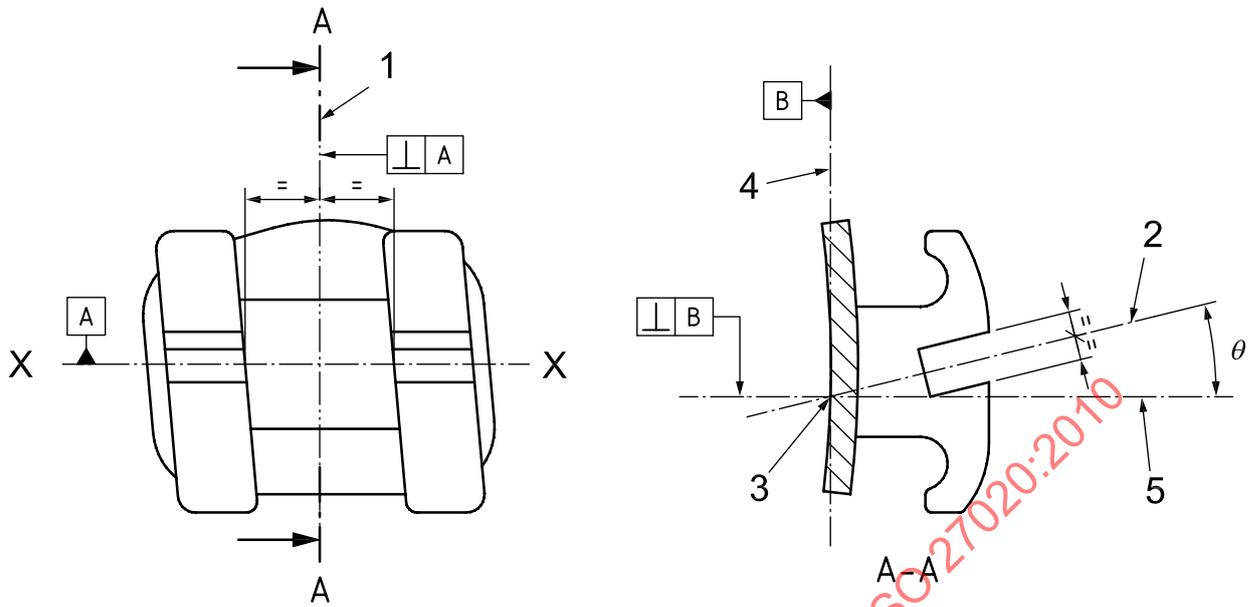

a) Angle of torque in-face

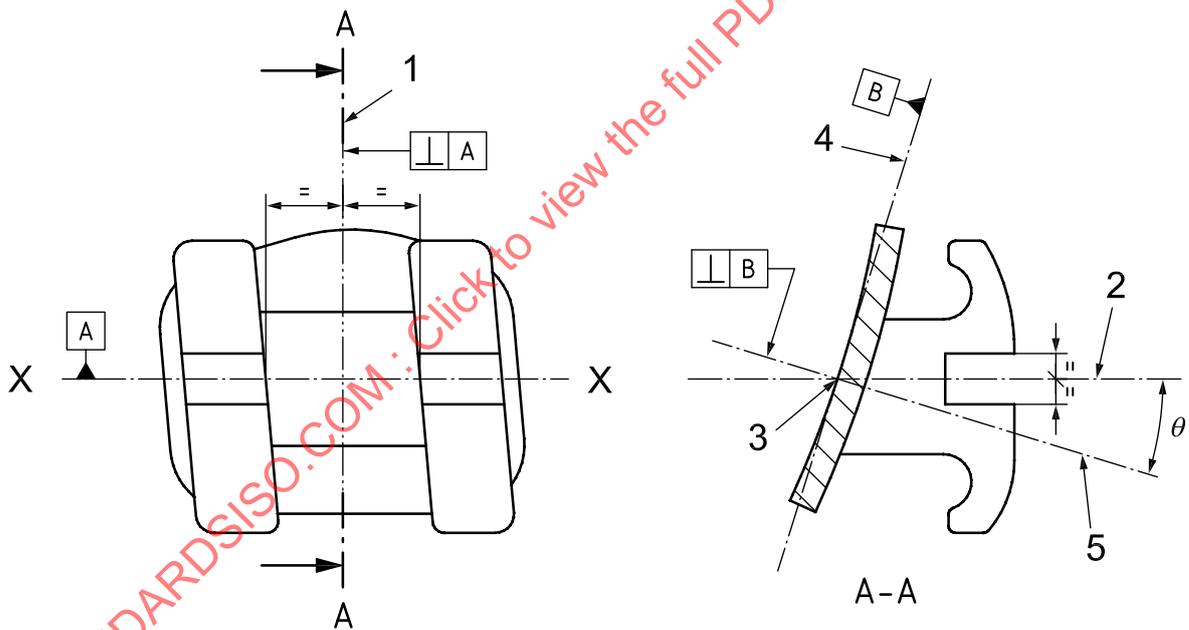

b) Angle of torque in-base

Key


- X mesial-distal long-axis of the slot
- 1 plane perpendicular to the X-axis in the centre of the slot
- 2 line bisecting the slot in the occlusal-gingival direction
- 3 point of intersection of line 2 to the tooth side surface of the base
- 4 line tangent to the tooth side surface of the base passing through point 3
- 5 line perpendicular to line 4 passing through point 3
- θ angle of torque

Figure 1 — Angle of torque of the bracket or tube

5.2.2.2 **Angulation** (see Figure 2)

Viewing from the labial/buccal, perform the following.

- a) Construct a line along the mesial-distal long-axis of the slot (X-axis).
- b) Construct line 1 (see A–A in Figure 2) perpendicular to the X-axis in the centre of the slot (depicted by the symbol $\perp A$).
- c) Construct line 2 along the central occlusal-gingival axis.
- d) Measure the angulation, α , between line 1 and line 2.

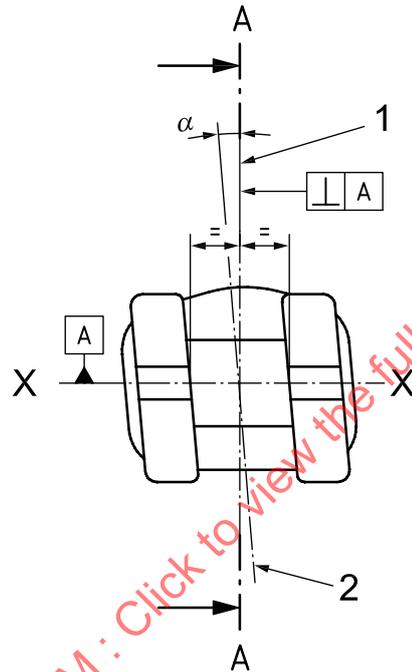

Key

- X mesial-distal long-axis of the slot
- 1 line perpendicular to the X-axis in the centre of the slot
- 2 line along the central occlusal-gingival axis
- α angulation

Figure 2 — Angulation of the bracket or tube

5.2.2.3 **In-out** (see Figure 3)

5.2.2.3.1 **Bracket in-out** [see Figure 3 a)]

When viewed from the labial/buccal, perform the following.


- a) Construct a line along the mesial-distal long-axis of the slot (X-axis).
- b) Construct plane 1 [see A–A in Figure 3 a)] perpendicular to the X-axis in the centre of the bracket slot (depicted by the symbol $\perp A$).

- c) When viewed on plane 1, construct plane 2 bisecting the slot in the occlusal-gingival direction.
- d) In plane 1, measure the in-out distance l_i from the bottom of the slot to the intersection of plane 2 with the tooth side surface of the base.

5.2.2.3.2 Tube in-out [see Figure 3 b)]

When viewed from the gingival, perform the following.

- a) Construct plane 2 bisecting the slot of the tube in the occlusal-gingival direction.
- b) When viewed on plane 2, construct plane 3 along the mesial end-face of the slot of the tube.
- c) Measure along the intersection of planes 2 and 3, then measure the in-out distance l_i from the bottom of the slot to the tooth side surface of the base.

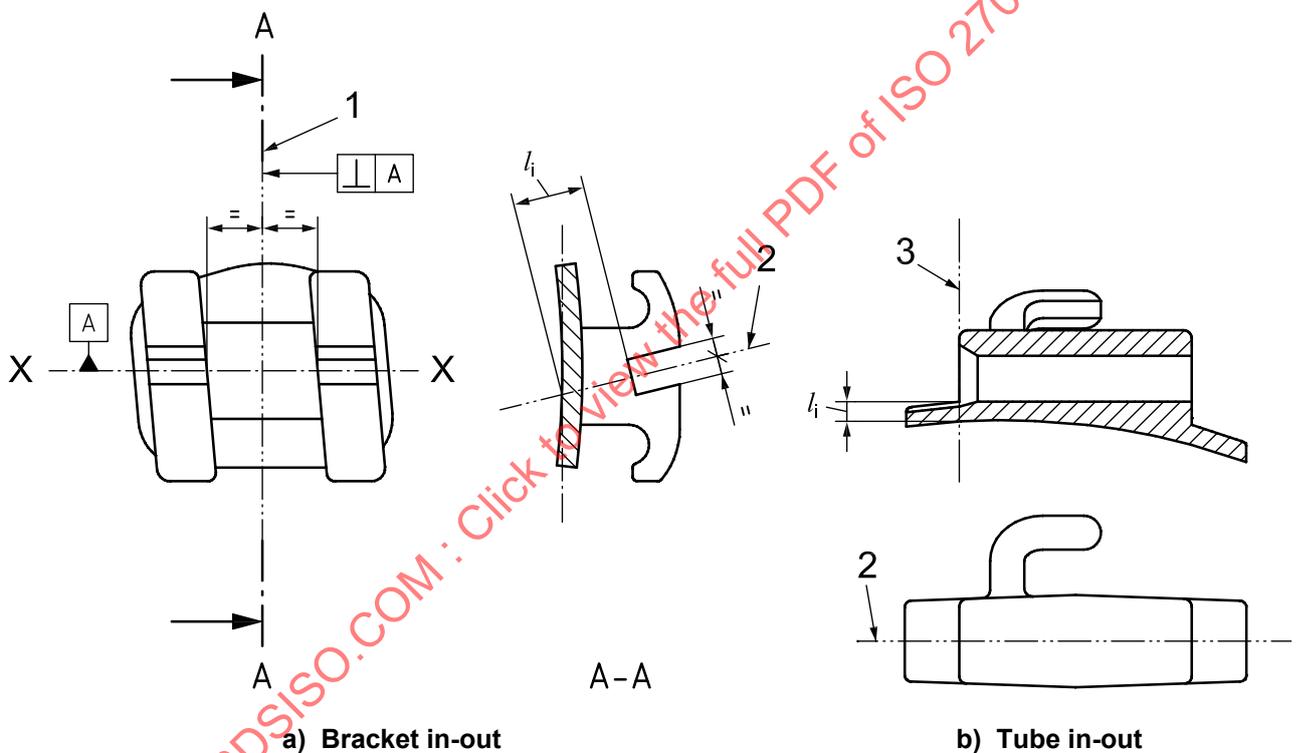

Key


- X mesial-distal long-axis of the slot
- 1 plane perpendicular to the X-axis in the centre of the bracket slot
- 2 plane bisecting the slot in the occlusal-gingival direction
- 3 plane perpendicular to an archwire that is fully engaged within the slot at the mesial end-face of the tube
- l_i in-out

Figure 3 — In-out

5.2.2.4 Rotational offset (see Figure 4)

When viewed from the occlusal, perform the following.

- a) Construct line 1 parallel to the direction of a straight section of wire that fills the slot.
- b) Construct lines 2 and 3 along the mesial and distal end-faces toward the base.
- c) Construct line 4 connecting the points of intersection of lines 2 and 3 with the tooth side surface of the base.
- d) Measure angle δ between lines 1 and 4 as rotational offset.

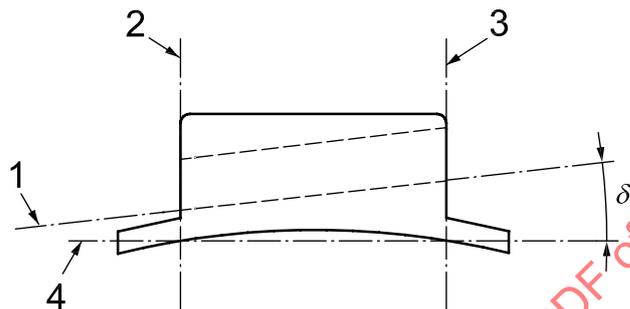

Key

- 1 line parallel to the floor of the slot
- 2, 3 lines along the mesial and distal end-faces of the slot, respectively
- 4 line connecting the points of intersection of lines 2 and 3 at the base
- δ rotational offset

Figure 4 — Rotational offset — Occlusal view

5.2.2.5 Slot depth (see Figure 5)

When viewed from the labial/buccal, perform the following.

- a) Construct a line along the mesial-distal long-axis of the slot (X-axis).
- b) Construct plane 1 (see A–A in Figure 5) perpendicular to the X-axis in the centre of the bracket slot (depicted by the symbol $\perp A$).
- c) When viewed on plane 1, project line 2 tangent to the shorter side of the intended slot and the buccal/labial side of the rectangle.
- d) Measure the slot depth as the dimension d .

5.2.2.6 Slot height (see Figure 5)

When viewed from the labial/buccal, perform the following.

- a) Construct a line along the mesial-distal long-axis of the slot (X-axis).
- b) Construct plane 1 (see A–A in Figure 5) perpendicular to the X-axis in the centre of the bracket slot (depicted by the symbol $\perp A$).
- c) When viewed on plane 1, measure the maximum occlusal-gingival dimension h .